


Personal floor scales KERN MPK · MPL

Personal scales with column for convenient, ergonomically-optimised operation




KERN MPK 200K-1P

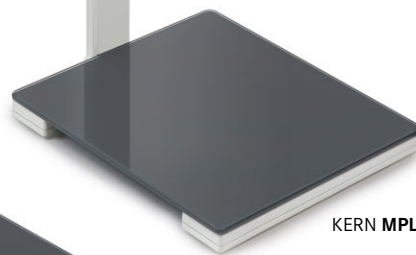

KERN MPL 200K-1P

Your KERN specialist dealer

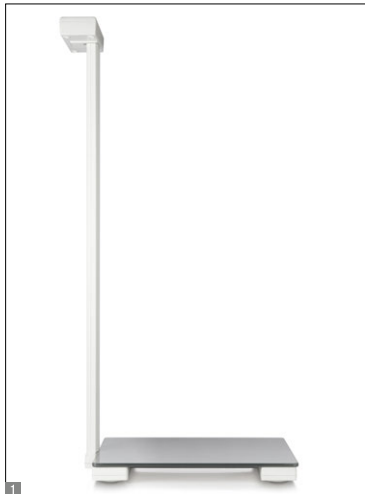

Features

- This personal scale is the ideal measuring instrument for elderly care, paramedics, fitness clubs or in the private sector and can handle weighing overweight persons weighing up to 250 kg
- The display device is easy to operate through the column which is firmly fixed to the scale and the weight can be read easily, even under poor lighting conditions
- Sturdy glass weighing plate for secure, comfortable standing position
- For easy and hygienic cleaning
- Particularly flat design facilitates access
- Easy to operate:
simply tap on it – stand on it – read the weight
- Integrated AUTO-OFF function to preserve the batteries
- KERN MPL: BMI function to determine underweight/normal weight/surplus weight
- Weight displayed in kg or lb
- Large high-contrast display that is easy to read
- KERN MPL: Adjusting function and tare function to zero the balance
- Delivered in appealing packaging
- Adjusting program CAL for quick setting of the balance accuracy using an external test weight

Technical data

- LCD display, digit height 28 mm
- Dimensions of weighing plate W×D×H 350×350×45 mm
- Overall dimensions (incl. stand) W×D×H 350×420×955 mm
- Batteries included, 4×1,5 V AA, operating time up to 20 h
- Net weight approx. 3,8 kg

Accessories

- Cleaning cloths, alcohol-free cloths for wipe disinfectant, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive medical products. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01

STANDARD



Model	Weighing range	Readability	
	[Max] kg	[d] kg	
KERN			
MPK 200K-1P	250	0,1	
MPL 200K-1P	250	0,1	

Pictograms



Adjusting program CAL:

For quick setting up of the balance's accuracy.
External adjusting weight required



Weighing units:

Can be switched to e.g. nonmetric units at the touch of a key. Please refer to website for more details



Battery operation:

Ready for battery operation. The battery type is specified for each device



Weighing principle: Strain gauges

Electrical resistor on an elastic deforming body



Package shipment:

The time required for internal shipping preparations is shown in days in the pictogram