

SC-240

PROFESSIONAL BODY COMPOSITION ANALYZER





The world's first portable, high capacity, professional grade, FDA Cleared body composition monitor.

The perfect tool to monitor athletes training for an event, or recovering from injury

Battery powered with convenient built-in carry handle

USB output is EMR ready, providing ready capability for HITECH Act compliance

Performs rapid analysis of clients up to 440 pounds

The most cost effective system providing professional results, from Tanita THE STANDARD

Expanded capability with Tanita Health Ware[®] Software, the best value



TANITA THE STANDARD

SC-240 **PROFESSIONAL GRADE MONITOR - WEIGHT AND BODY COMPOSITION SCALE**

This FDA cleared body composition monitor is light and highly portable with a built-in handle, yet fully featured for all age ranges with a 440 lb weight capacity. Weighing just 10 lb, the large platform measures weight and instantly calculates body fat% and total body water% in addition to BMI with clinically proven accuracy onto a large LCD display. Offering in-depth individual analysis of body composition with clinically proven accuracy. The monitor also has a USB port allowing the highly repeatable data to be easily transferred to Tanita Health Ware® Software (sold separately). Tanita Health Ware® Software provides custom data recording, tracking and analysis for an unlimited number of users. Reports are easily created using standard color printers, no special forms are required. The SC-240 and Tanita Health Ware® combination provides detailed analysis and long term trend monitoring of body composition data, vital in assessing the effectiveness of any health or fitness program.

Features:

Data Shown on the LCD Display:

- \checkmark Weight
- Body Fat %
- Body Water % \checkmark
- \checkmark BMI

Additional data available using Tanita Health Ware®:


- Fat mass \checkmark
- \checkmark Fat free mass
- Body water mass \checkmark
- Muscle mass \checkmark
- \checkmark Bone mass
- \checkmark Visceral fat rating
- \checkmark Metabolic age
- \checkmark Muscle score
- \checkmark Daily caloric intake (DCI)
- Basal metabolic rate (BMR) \checkmark
- Ideal body weight \checkmark
- \checkmark Degree of obesity
- \checkmark BMR score
- \checkmark Impedance

Specifications:

- \checkmark Capacity: 440lb (200 kg)
- Weight Graduation: 0.2lb (100 g) \checkmark
- ✓ Body Fat Increments: 0.1%
- ✓ USB Type B Output
- \checkmark Dimensions: 17.2" x 13.4" x 2.1"
- Weight: 10.0lbs \checkmark
- ✓ Power Supply: AC Adapter (included)
 - or 6 x AA (not included)
- 3 vear warrantv



Tanita HealthWare™ detailed Body Composition Report

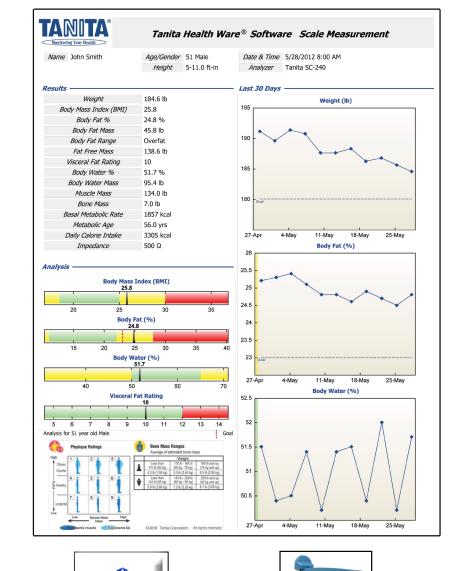

2625 S. Clearbrook Drive / Arlington Heights, IL 60005 / Phone: 847-640-9241 / Fax: 847-640-7978 / E-mail: 4health@tanita.com / Web: www.Tanita.com Tanita Health Ware® is a registered trademark of TANITA Corporation in the United States

Optional

Accessories:

C-300

Carrying Case

HR-200

Height Rod